

Customer No.	026418	
IN THE UNITED STATES PATENT AND TRADEMARK OFFICE		
Attorney's Docket No.:	GK-ZEI-3254 / 500343.20274	
U.S. Application No.:	10/511,483	
International Application No.:	PCT/EP03/04024	
International Filing Date:	APRIL 17, 2003	17 APRIL 2003
Priority Date Claimed:	APRIL 17, 2002	17 APRIL 2002
Title of Invention:	METHOD FOR THE SPECTROSCOPIC DETERMINATION OF THE OXYGEN SATURATION OF BLOOD IN THE PRESENCE OF OPTICAL DISTURBANCE VARIABLES	
Applicant(s) for (DO/EO/US):	Dietrich SCHWEITZER, Martin HAMMER and Eike THAMM	

**VERIFIED
ENGLISH
TRANSLATION
OF
PUBLICATION NO.
WO 03/086193**

Published on October 23, 2003

VERIFICATION OF A TRANSLATION

I, the below named translator, hereby declare that:

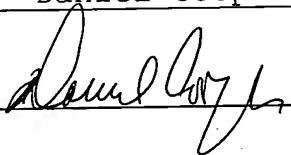
My name and post office address are as stated below:

That I am knowledgeable in the English language and in the language in which the below-identified international application was filed, and that I believe the English translation of the international application PCT/EP03/04024 is a true and complete translation of the above-identified international application as filed.

I hereby declare that all statements made herein of my own knowledge are true and that all statements made on information and belief are believed to be true; and further that these statements were made with the knowledge that willful false statements and the like so made are punishable by fine or imprisonment, or both, under section 1001 of Title 18 of the United States Code and that such willful false statements may jeopardize the validity of the application or any patent issued thereon.

Date December 16, 2004

Full name of the translator Daniel Cooper

Signature of the translator 

Post Office Address 1310 Felicity Street, New Orleans, LA 70130

10/511483

Docket No.: GK-ZEI-3254/500343.20274

METHOD FOR THE SPECTROSCOPIC DETERMINATION OF
THE OXYGEN SATURATION OF BLOOD IN THE PRESENCE OF OPTICAL
DISTURBANCE VARIABLES

[0001] The invention is directed to a method for the spectrometric determination of the oxygen saturation of blood in the presence of optical disturbance variables such as those also presented by pigmented and light-scattering surrounding tissue and/or also the vascular wall itself. This problem of determining the oxygen saturation of blood while excluding as far as possible the influence of these factors affecting measurement accuracy occurs particularly in noninvasive, in-vivo or in-vitro examinations of blood vessels located in front of, behind, or in said pigmented and scattering tissue, for example, in the examination of blood vessels of the ocular fundus or other tissue areas in the body such as the skin and organs that are accessible by endoscopy.

[0002] It is generally known that the absorption spectrum of the red blood pigment hemoglobin changes with oxygen saturation (e.g., van Assendelft, O. W., Spectrophotometry of hemoglobin derivatives, Assen: Royal Vangorcum, 1970). Accordingly, it is possible to determine the oxygen saturation of a hemoglobin sample by comparing the spectrum of the sample to the spectra of completely oxygenated and completely reduced hemoglobin.

[0003] Recent work in oxymetry in the ocular fundus using the Lambert-Beerschen law, i.e., taking only absorption into account, has been published by Smith et al. and others (Smith, M. H., Denninghoff, K. R., Lompado, A., Hillman, L. W., Effect of multiple light paths in retinal vessel oxymetry, *Appl. Opt.* **39**, 2000, 1183-1193). Numerous patented arrangements and methods are based on this principle (e.g., US 4,485,820; US 5,119,814; US 5,308,919; US 4,253,744; US 4,305,398; US 5,776,060; US 5,935,076; DE 199 20 157 A1; US 5,318,022).

[0004] However, the hemoglobin does not exist in isolation in in-vivo measurement, but is enclosed in the erythrocytes. The scattering of light on the erythrocytes has a considerable influence on the extinction spectrum of the blood. However, based on the findings of the multiple scattering theory of Twersky (Twersky, V., Absorption and multiple scattering by

biological suspensions, *J. Opt. Soc. Amer.* **60**, 1970, 1084-1093), the influences of scattering and absorption can be separated. On this basis, Pittman and Duling describe a method for determining the oxygen saturation in whole blood from measurements taken in transmission at a wavelength of 555 nm and at isosbestic points at 522 nm and 546 nm (Pittman, R. N., Duling, B. R., A new method for the measurement of percent oxyhemoglobin, *J. Appl. Physiol.*, **38**, 1975, 315-320). This method has been used by Delori to determine oxygen saturation in retinal vessels (Delori, F. C., Noninvasive technique of oximetry of blood in retinal vessels, *Appl. Opt.* **27**, 1988, 113-1125).

[0005] However, investigations by Hammer et al. (Hammer, M., Leistritz, S., Leistritz, L., Schweitzer, D., Light paths in retinal vessel oxymetry, *IEEE Trans Biomed Eng* **48** (5), 2001, 592-8) show that the reflection spectra measured on retinal vessels are influenced not only by the absorption of the hemoglobin and the scattering in the blood and in the tissue surrounding the vessels, but also by the melanin located in the retinal pigment epithelium and in the choroid. This is also true of vessels in the skin or other pigmented organs.

[0006] Correcting falsification of the hemoglobin spectra by means of other chromophores and correction for spectroscopic oxymetry have been attempted in previous literature by scaling the spectra measured on a vessel to measurements next to the vessel (e.g., DE 199 20 157 A1; US 5,935,076; Delori, F. C., Noninvasive technique for oximetry of blood in retinal vessels, *Appl. Opt.* **27**, 1988, 113-1125; Schweitzer, D., Hammer, M., Kraft, J., Thamm, E., Königsdörffer, E., Strobel, J., In Vivo Measurement of the Oxygen Saturation at the Normal Human Eye, *IEEE Trans. Biomed. Eng.* **46**, 1999, 1454-1465). However, this approach does not take into account the extremely complicated relationships (Hammer, M., Leistritz, S., Leistritz, L., Schweitzer, D., Light paths in retinal vessel oxymetry, *IEEE Trans Biomed Eng* **48** (5), 2001, 592-8) of the beam propagation in the blood vessel and the tissue surrounding it.

[0007] The exact light propagation in the biological tissue cannot always be fully described physically. Even efforts to recreate these processes for eliminating disturbances in the most comprehensive and realistic manner possible (DE 199 20 157 A1) and to model the optics of the living or dead biological tissue surrounding the blood vessel (DE 44 33 827 A1) have not led to more exact measurements than the methods mentioned above, which are already relatively time-consuming and require intensive computation. In view of this, the

methods are only conditionally suitable especially for routine examinations and screening examinations. In particular, the determination of oxygen saturation at every point of a two-dimensional, graphic recording, which is important for clinical practice, requires a method that is fast on the one hand and that compensates for optical and spectrometric disturbances due to the tissue environment on the other hand.

[0008] Therefore, it is the object of the invention to determine the oxygen saturation with high accuracy by a method which is as simple and as fast as possible.

[0009] In order to meet the above-stated object, spectral measurements are generated by transmission measurement and reflection measurement in a measurement spectrum at wavelengths that are isosbestic for hemoglobin and oxyhemoglobin and at least one other measurement value is generated at a wavelength at which the reference values of hemoglobin and oxyhemoglobin differ, and these measurements are compared with known reference values of the reference spectra of hemoglobin and oxyhemoglobin in that:

[0010] a) at least two said spectral measurement values (M_{i1} , M_{i2}) at wavelengths (λ_{i1} , λ_{i2}) that are isosbestic for hemoglobin and oxyhemoglobin and at least the other measurement value (M_a) at a wavelength (λ_a) at which the reference values of hemoglobin and oxyhemoglobin differ as far as possible in the reference spectra are detected in the measurement spectrum, wherein an auxiliary function (F_H) is generated at least from two of the measurement values (M_{i1} , M_{i2}) for isosbestic wavelengths (λ_{i1} , λ_{i2}),

[0011] b) a reference function (F_R) is generated in the reference spectra from the reference values (R_{i1} , R_{i2}) corresponding to the at least two measurement values (M_{i1} , M_{i2}) determined in the measurement spectrum for the same isosbestic wavelengths (λ_{i1} , λ_{i2}) of hemoglobin and oxyhemoglobin, which reference function (F_R) is of the same type,

[0012] c) a correction function (F_K) is generated from the auxiliary function (F_H) of the measurement spectrum in which said at least two measurement values (M_{i1} , M_{i2}) lie for isosbestic wavelengths (λ_{i1} , λ_{i2}) and from the reference function (F_R) of the reference spectra in which the at least two reference values (R_{i1} , R_{i2}) corresponding to the at least two measurement values (M_{i1} , M_{i2}) lie, and a corrected auxiliary function (F_{HK}) identical to the reference function (F_R) in the reference spectra is generated in a corrected measurement spectrum by means of this correction function (F_K), and

[0013] d) the oxygen saturation of the blood is determined from the other measurement value (M_a'') converted to the corrected auxiliary function (F_{HK}) of the corrected measurement spectrum in relation to the reference values for hemoglobin and oxyhemoglobin at this wavelength (λ_a).

[0014] It is advantageous when the spectral measurement values and reference data are determined logarithmically and the auxiliary function and reference function are formed by a straight line on which two of the measurement values or reference values lie for isosbestic wavelengths. The correction function formed from the linear auxiliary function and reference function can also give a linear corrected auxiliary function of the corrected measurement spectrum. A constant multiplier is applied to the rest of the spectral measurement values, i.e., the spectral measurement value for the third isosbestic wavelength and the other measurement value at a wavelength at which the reference values of hemoglobin and oxyhemoglobin differ as far as possible in the reference spectra. This constant multiplier is determined in such a way that the third isosbestic measurement value of the corrected measurement spectrum that is corrected by scaling conforms to the corresponding reference value of the reference spectra. In this specific case, the difference between the reference values for hemoglobin and oxyhemoglobin can be scaled linearly between 0 and 1. The oxygen saturation of the blood is determined in relation to this scale from the other measurement value that is converted to the corrected auxiliary function of the corrected measurement spectrum.

[0015] For purposes of a clear two-dimensional representation of the oxygen saturation of the blood, four monochromatic individual images are generated at said wavelengths and the oxygen saturation is determined for each image point.

[0016] Surprisingly, it has been shown that in comparison with the examination methods mentioned in the introduction, the oxygen saturation of the blood can be determined by the proposed method with the same accuracy but with substantially less effort (minimally, four measurement values are required) for measurement, calculation and evaluation. The method makes possible a two-dimensional, spatially dependent representation of the measurements allowing a manageable and fast evaluation. Only a few measurement values are required and only linear transformations are used. With these advantages of low expenditure and time-

efficient measurement evaluation, the proposed method is also suitable for screening examinations as well as for routine and early detection examinations.

[0017] The invention will be described more fully in the following with reference to an embodiment example shown in the drawing.

[0018] Fig. 1 is a graph showing measurement values and reference values in logarithmic form in the wavelength range between 400 nm and 700 nm; including three isosbestic wavelengths $\lambda_{i1} = 522$ nm, $\lambda_{i2} = 586$ nm, $\lambda_{i3} = 569$ nm and the other wavelength $\lambda_a = 555$ nm;

[0019] Fig. 2 is a graph with spectral reference values according to Fig. 1 showing linear reference function F_R and corrected measurement values M' ;

[0020] Fig. 3 is a graph with spectral reference values according to Fig. 1 showing linear reference function F_R with corrected measurement values M'' and with scaling for reading off oxygen saturation.

[0021] The reflection or transmission of tissue at a point or with spatial resolution in an image is measured at three isosbestic wavelengths λ_{i1} , λ_{i2} and λ_{i3} ($\lambda_{i1} = 522$ nm, $\lambda_{i2} = 586$ nm, $\lambda_{i3} = 569$ nm) as measurement data M_{i1} , M_{i2} and M_{i3} and at a different wavelength λ_a (555 nm), at which the absorption coefficients of oxygenated and reduced hemoglobin differ, as measurement value M_a and is compared with the reflection or transmission of hemoglobin or whole blood with oxygen saturations of 0% and 100%, respectively, at these wavelengths (reference data R_1 , R_2 , R_3 , $R_a^{0\%}$ and $R_a^{100\%}$) according to the following method:

[0022] 1. All measurement data and reference data are logarithmized. Fig. 1 shows the measurement data and reference data represented logarithmically, which are, in this example, the reflection of a retinal vein (measurement data) and the transmission of a layer of whole blood with a thickness of 0.1 mm (reference data). For the sake of clarity, the complete spectra between 400 nm and 700 nm are shown. The wavelengths $\lambda_{i1} = 522$ nm, $\lambda_{i2} = 586$ nm, $\lambda_{i3} = 569$ nm and $\lambda_a = 555$ nm used in this example are shown.

[0023] 2. A linear auxiliary function F_H of the wavelength is calculated in the measurement spectra such that its values at the isosbestic wavelengths λ_{i1} and λ_{i2} agree with the measurement data M_{i1} and M_{i2} at these wavelengths.

[0024] 3. A linear reference function F_R of the wavelength is calculated in the reference spectra in such a way that its values at the isosbestic wavelengths λ_{i1} and λ_{i2} agree with the reference data R_{i1} , R_{i2} at these wavelengths.

[0025] 4. The measurement data are corrected additively at each wavelength by the difference of the linear functions F_H and F_R in such a way that at the isosbestic wavelengths λ_{i1} and λ_{i2} they agree with the reference data: $M'_\lambda = M_\lambda + F_R - F_H$.

[0026] Fig. 2 again shows the reference data, the linear reference function F_R of the wavelength and the corrected measurement data M' . This correction compensates for extinctions which exist in addition to the absorption of the hemoglobin and whose spectra in the wavelength range of 522 nm to 586 nm can be assumed or approximated as linear in the logarithmic scale. The embodiment example under consideration concerns absorptions of melanin and the anterior optical media and the scattering in the blood and in the surrounding tissue.

[0027] 5. The corrected measurement data M' are scaled by a factor of less than or greater than 1 around the linear reference function F_R in such a way that they agree with the reference value R_{i3} at the isosbestic wavelength λ_{i3} :

$$M''_\lambda = F_R(\lambda) + \frac{(M'_\lambda - F_R(\lambda))(R_{i3} - F_R(\lambda_{i3}))}{M'_{i3} - F_R(\lambda_{i3})}$$

[0028] Fig. 3 shows the spread (or scaling) of the corrected measurement data M' resulting in M'' around the linear reference function F_R which is carried out in such a way that corrected measurement data and reference data agree at the isosbestic wavelength λ_{i3} (569 nm). This correction compensates for different absolute values of the measurement data and reference data that occur due to different illumination conditions and measurement conditions.

6. The position of M''_a on a scale contained linearly between

$R_a^{0\%}$ and $R_a^{100\%}$ indicates the oxygen saturation OS :

$$OS = \frac{M''_a - R_a^{0\%}}{R_a^{100\%} - R_a^{0\%}}$$

[0029] The scaled readout of the oxygen saturation (OS) between the values 0 and 1 is likewise shown in Fig. 3. The oxygen saturation OS in the embodiment example is 0.69.

Reference Numbers

M	-	measurement data
M'	-	measurement data corrected by the addition of the corrected auxiliary function
M''	-	measurement data M' corrected by applying a factor
F _R	-	reference function
OS	-	oxygen saturation
λ	-	wavelength